FAZZINI

electrosurgery

7

Elettrochirurgia
Electrosurgery
Bistouris electriques et acc.
Electrocirugia
Elektrochirurgie

ELETTROBISTURI
ELECTROSURGICAL UNITS
BISTOURIS ELECTRIQUE
ELECTROBISTURI
ELEKTROCHIRURGICKÉ JEDNOTKY
50/80/120/160/200/240/380 W

PINZE BIPOLARI
BIPOLAR FORCEPS
PINCES BIPOLAIRES
PINZAS BIPOLARE
BIPOLÁRNÍ PINZETY

ELETTRODI
ELECTRODES
ELECTRODES
ELECTRODOS
ELEKTRODY

MANIPOLI
ELECTROSURGERY PENCILS
MANCHES PORTE-ELECTRODES
MANGOS PORTA ELECTRODO
ELEKTROCHIRURGICKÉ DRŽÁKY NÁSTROJŮ

PIASTRE ~ CAVI DI CONNESSIONE
GROUND ELECTRODE ~ CONNECTION CABLE
ELECTRODES NEUTRES ~ CABLES
ELECTRODO NEUTRO ~ CABLE DE CONEXIÓN
NEUTÁLNÍ ELEKTRODY, PROPOJOVACÍ KABELY

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DIATROM® 50D ELECTROSURGICAL UNIT
FOR MONOPOLAR LIGHT SURGERY

DIATROM® 50D is suitable to light monopolar pure cut as well as cut/coagulation and coagulation.
Under microcontroller watch it has the maximum output linearity, safety control allowing long lasting trouble free operation.
The operative function controlled by keys and the digital display of the output power makes the surgery's work simple and reliable.

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DIATROM® ELECTROSURGICAL UNIT
MONO/BIPOLAR LIGHT SURGERY

DIATROM® is a radiosurgical unit which is suitable to microcoagulative operations as well as to light and medium monopolar and bipolar coagulators. DIATROM® by means of the selection of the unit's performance, allows CUT, BLEND, COAG FORCED, COAG SOFT, BIPOLAR with special adaptor .
The digital reading of the supplied power makes the applications safe and highly reliable. DIATROM® allows a highly professional surgery thanks to the ergonomic and safety solutions normally utilised in equipment with higher output power used in the surgical room.
In fact the neutral electrode connection is continuously moni-tored and the presence of the controls of the performance and of the delivery of the
power placed on the handpiece, they allows the making of the operations without turning away the surgeon attention from the operating theatre.

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DIATROM® 200

DIATROM 200® is a electrosurgical unit suitable to effect middle monopolar or bipolar electrosurgery. **DIATROM® 200** by means of the selection of the unit's performance, allows CUT, CUT 90%, CUT 80%, CUT 70%, BLEND, COAG FORCED, COAG SOFT, BIPOLAR CUT, BIPOLAR COAG, with special adaptor. The planning of the power output and pause times allow the use for collecting bioptic samples, with reduced modification of the withdrawn sample and reduced thermal damage to the tissue object of the collection. Besides, it allows stopping the coagulation with monopolar or bipolar forceps automati-cally. **DIATROM® 200** allows the accurate predisposition of the output power level, with resolution of 1 Watts and due to programmable time of supply it is usable for interventions of micro-surgery and precise dermatologic surgery like treatment of telangiectases, spider angiomas, vascular nevi, etc. The coagulation can be performed with high crest factor cur-rent to get an elevated thermal effect or, with low crest factor current, to get coagulation in absence of carbonisation. The control of the unit is made through pushbuttons and indicators set on the frontal panel.
The type of surgical operations that may be performed are all those that require the monopolar or bipolar cut or coagula-tion. When bipolar forceps have to be used the special adapter (REF 00498.04) is required. 9 different programs of work related to 9 different opera-tional needs could be stored and recalled simply from the operator.

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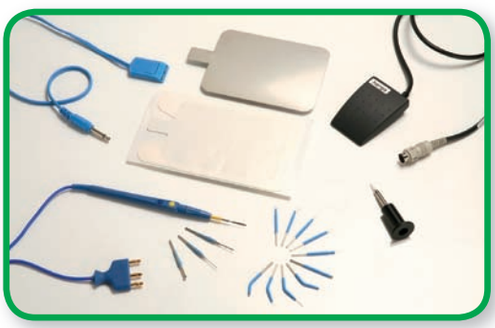
DIATROM® FLASH ~ PROGRAMMABLE ELECTROSURGICAL UNIT
MONO/BIPOLAR

DIATROM® FLASH is a diathermo-coagulator unit which is suitable to microsurgery as well as to precision monopolar and bipolar surgery.
DIATROM® FLASH can supply RF operative power either for programmable time as short as thousends of a second or for the whole driving time carried out by a handpiece or a footswitch.
DIATROM® FLASH by means of the selection of the unit's performance, allows CUT, BLEND, COAG FORCED, COAG SOFT, BIPOLAR with special adaptor.
DIATROM® FLASH allows a highly professional surgery thanks to the ergonomic and safety solutions nor-mally utilised in equipments with higher output power used in the sur-gical room.

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DIATROM® FLASH ~ HF 200 W

DIATROM® FLASH HF 200 W is a HF electrosurgical unit suitable to effect middle monopolar or bipolar electrosurgery. The equip-ment can be used with power emission time equal to that of actuation of the output circuit as well as with emission time preset from 1ms up to 999 ms. The programmed pulse can be single or repeated from 1 to 250 times per second. **DIATROM® FLASH HF 200 W** by means of the selection of the unit's performance, allows CUT, CUT 90%, CUT 80%, CUT 70%, BLEND, COAG FORCED, COAG SOFT, BIPOLAR, with special adaptor. The planning of the power output and pause times allow the use for collecting bioptic samples, with reduced modification of the withdrawn sample and reduced thermal damage to the tissue object of the collection. Besides, it allows stopping the coagulation with monopolar or bipolar forceps automatically. **DIATROM® FLASH HF 200 W** allows the accurate predisposition of the output power level, with resolution of 1 Watts and due to programmable time of supply it is usable for interventions of micro-surgery and precise dermatologic surgery like treatment of telangiectases, spider angiomas, vascular nevi, etc. The coagulation can be performed with high crest factor cur-rent to get an elevated thermal effect or, with low crest fac-tor current, to get coagulation in absence of carbonisation. The control of the unit is made through pushbuttons and indicators set on the frontal panel.
The type of surgical operations that may be performed are all those that require the monopolar or bipolar cut or coagula-tion. When bipolar forceps have to be used the special adapter (REF 00498.03) is required. 9 different programs of work related to 9 different opera-tional needs could be stored and recalled simply from the operator.

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DIATROM® 240
HIGH FREQUENCY ELECTROSURGERY EQUIPMENT
240 WATT

EASY TO USE: The operation modality and the output controls selection is made simply by means of keys provided with iconographic description. All the indications of operative condition and of selected output level are digital bright.
CUT: The monopolar cut can be performed by pure sine current to obtain an incision with a minimum damage for the adjacent tissues and sure recovery for first intention. Two different modality of modulation allows to reduce the bleeding during the monopolar cut. In the bipolar modality the cut can be effected with pure sine current.
HEMOSTASIS: When the surgical action requires to obliterate a bloody vessel the availability of current with different coagulative power allows to effect this type of work both with monopolar and bipolar surgery.
FULGURATION: When it is necessary to cause a superficial necrosis of wide areas, protecting the deeper tissues, it is available for the operation a special coagulant current that dehydrates the superficial layer of the tissue and reduces the possibility of damage those underlying.
OUTPUT CONTROL: Output power levels are adjustable by means of incremental encoders that allow the exact predetermination and reproductivity of the values. The output function is linearized and managed by means of ultramodern electronic microcontrollers which prevent the delivery when any condition of checked alarm happen.
SAFETY OF USE: Circuits dedicated to the safety of operator and patient had been particulary mind. That allows to maintain the leakage current well under the values pointed out from the international requirements, under every output condition. Monopolar and bipolar output circuits are completely floating and the use of the defibrillator is allowed while the reference plate patient is connected.
RELIABILITY: DIATROM ®240 uses all solid state electronic solid state circuits and it is realized with modular interchangeable unity for a fast and effective maintenance. Particular automatic control circuits protect the unities from damage also due to repeated improper use.
WORKABILITY: DIATROM® 240 can be used to execute monopolar and bipolar surgery in the many fields of application where it is requested high precision and reliability like in the endoscopic surgery.

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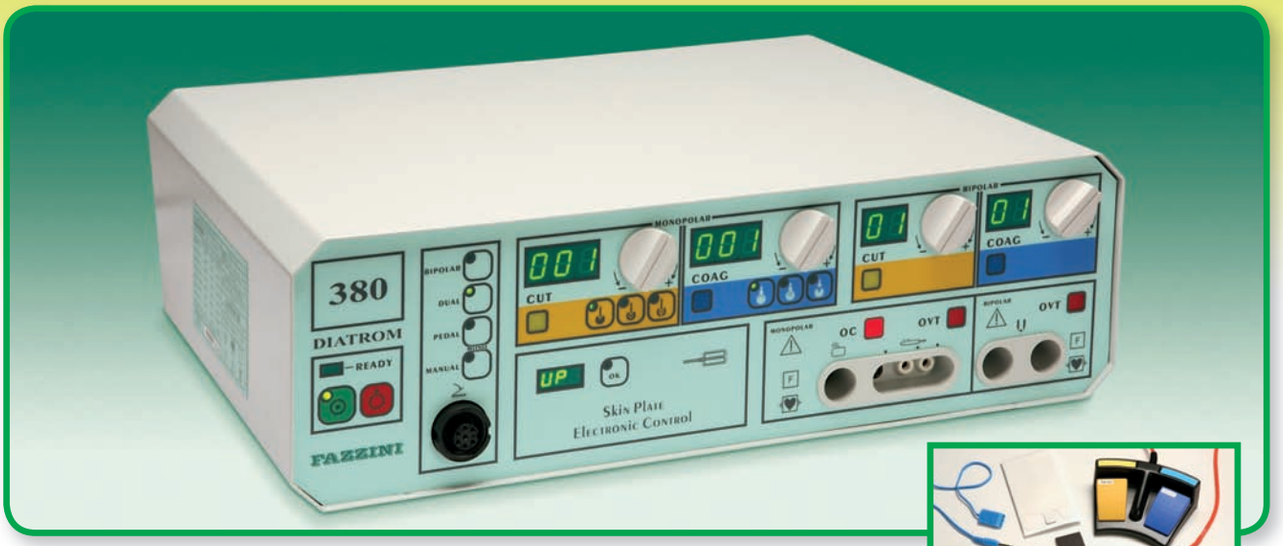
DIATROM® 380
HIGH FREQUENCY ELECTROSURGERY EQUIPMENT
380 WATT

EASY TO USE: The operation modality and the output controls selection is made simply by means of keys provided with iconographic description. All the indications of operative condition and of selected output level are digital bright.
CUT: The monopolar cut can be performed by pure sine current to obtain an incision with a minimum damage for the adjacent tissues and sure recovery for first intention. Two different modality of modulation allows to reduce the bleeding during the monopolar cut. In the bipolar modality the cut can be effected with pure sine current.
HEMOSTASIS: When the surgical action requires to obliterate a bloody vessel the availability of current with different coagulative power allows to effect this type of work both with monopolar and bipolar surgery.
FULGURATION: When it is necessary to cause a superficial necrosis of wide areas, protecting the deeper tissues, it is available for the operation a special coagulant current that dehydrates the superficial layer of the tissue and reduces the possibility of damage those underlying.
OUTPUT CONTROL: Output power levels are adjustable by means of incremental encoders that allow the exact predetermination and reproductivity of the values. The output function is linearized and managed by means of ultramodern electronic microcontrollers which prevent the delivery when any condition of checked alarm happen.
SAFETY OF USE: Circuits dedicated to the safety of operator and patient had been particulary mind. That allows to maintain the leakage current well under the values pointed out from the international requirements, under every output condition. Monopolar and bipolar output circuits are completely floating and the use of the defibrillator is allowed while the reference plate patient is connected.
RELIABILITY: DIATROM® 380 uses all solid state electronic solid state circuits and it is realized with modular interchangeable unity for a fast and effective maintenance. Particular automatic control circuits protect the unities from damage also due to repeated improper use.
WORKABILITY: DIATROM® 380 can be used to execute monopolar and bipolar surgery in the many fields of application where it is requested high precision and reliability like in the endoscopic surgery.

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TWIN OPERATIVE ACTIVITY
• Two output channels for monopolar cutting
• Two independent setting and output channels for the most frequent monopolar coagulation
CONTROLLED CUTTING PROCEDURES
• Enhanced cutting current for eschar free incision
• Cutting current with selectable coagulating property for bleeding controlled dissection
• Low capacitive coupling by Slow Blend current for endoscopic procedures
TWIN CONTROLLED COAGULATIVE PROCEDURES
• Speedy* coagulation and Deep**coagulation availableto two operators
• Spray coagulation available by hand switch or foot switch
WIDE BIPOLAR FEATURES
• Bipolar cutting current with selectable bleeding control
• Slow Blend cutting current
• Separately activable Automatic Start and Automatic Stop coagulation
SURGEON AND PATIENT TOTAL SAFETY
• Acceptability of contact impedance between patient body and return electrode continuously watched by SkinPlate Electronic Control (SPEC)
• Excessive power delivery time warned
• H.F.Leakage currents continuously monitored
• Mains voltage variations automatically compensated before to start the delivery of HF power
PERSONALIZED PROCEDURES
• Independent choise of coagulative currents and power levels by each of two operators
• Ten personal configurable and immediately recallableprocedures
EXTENTION DRIVER AND INTERLOCK SIGNALS
• Drived output connection to external Argon Plasma unitand Smoke Evacuator unit
AUTOMATIC SELF DIAGNOSIS
• Automatic diagnostic control of the actual complete functionality of the electronic unit and the connected accessories and coded report of found failures

** Speedy coagulation also said Fulgurate or Forced coagulation*
*** Deep coagulation als said Pin Point,Dessicate or Soft coagulation*

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Adjustment of the levels of power of the functions of the monopolar section CUT:
CUT, for pure cut or with coagulative properties, suited for dissection;
ENHANCED, for incision with reduced production of eschar;
BLEND, for endoscopic procedures minimally invasive,with low capacity coupling through the surrounding tissues.
Selection of the function to deliver.
Visualization of the level of power of the selected function,that can be activated through command of the connected accessories.

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Adjustament of the levels of power of the functions of the monopolar section COAG1:
SPEEDY, for rapid coagulation;
DEEP, for coagulation and deep dessication.
Selection of the function to deliver.
Visualization of the level of power of the selected function,that can be activated through command of the connected accessory.

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Observation of the state of mains power supply and electronic interrupters for lighting and turning off of the unit.The fluctuations of the main power supply voltage are automatically compensated before to start the delivery of the HF current.

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User interface for:
- Planning of the time of delay of the beginning of the automatic bipolar coagulation. This time can be programmed to 0,5s,1s or 1,5s.
- Reading of last 101 errors eventually happened.
- Recording of 10 personalized procedures; selection and recall of the same. In the personalized procedures it is possible to record both the selected monopolar and bipolar functions and the relative levels of power.

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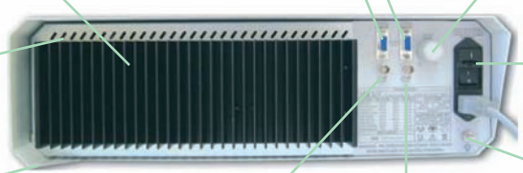
Adjustament of the levels of power of the functions of the monopolar section COAG2:
SPEEDY, for rapid coagulation;
DEEP, for coagulation and deep dessication;
SPRAY, for wide superficiale coagulation without contact;
ARGON SPRAY, supported by unity Argon-Plasma.
Selection of the function to deliver.
Visualization of the level of power of the selected function, that can be activated through command of the connected accessories.

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
Adjustament of the levels of power ofthe bipolar CUT:
CUT, for pure cut or with coagulative property, suited for dissection;
BLEND, low voltage current,proper to endoscopic operations, thanks to the low risk of power transfer by capacitive coupling to the surrounding tissues.
Digital visualization of the bipolar cut output power.

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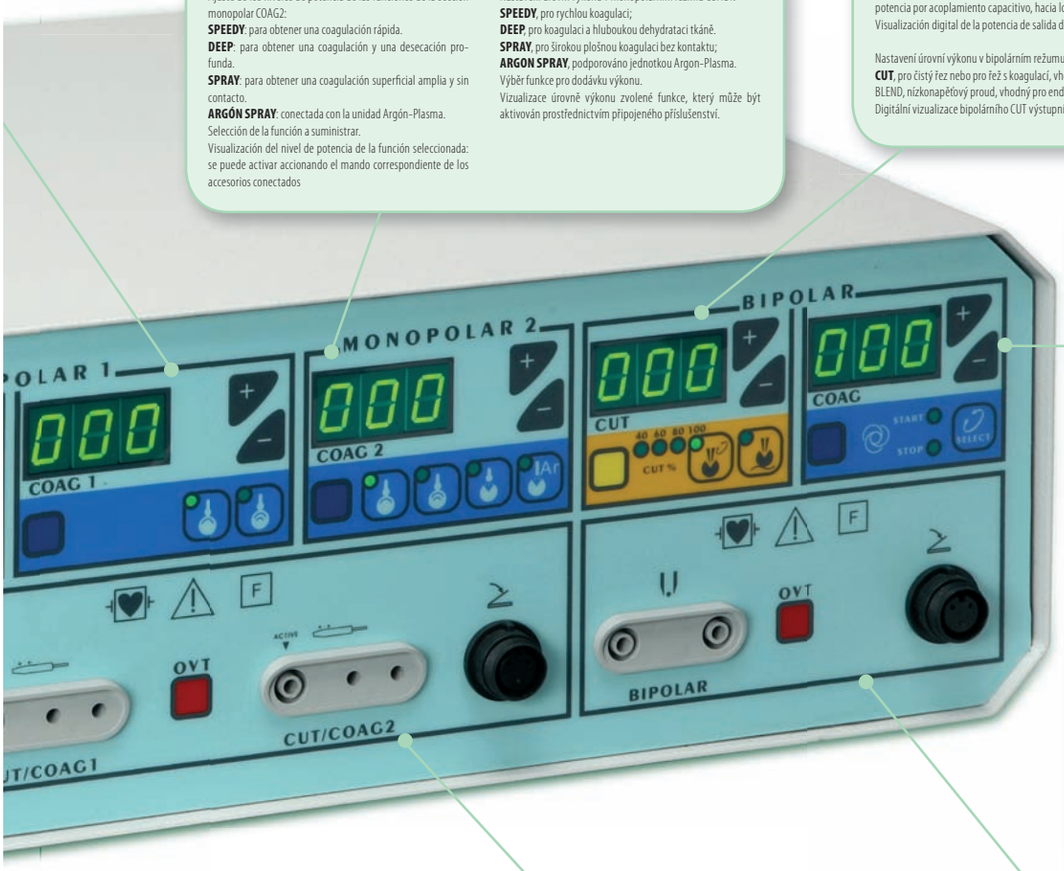Adjustment of the power level of bipolar coagulation:
Choise of the modality of delivery of the current can be also selected with automatic beginning (START) and/or automatic end (STOP). The automatic beginning of the delivery of current can take place0,5s, 1s or 1,5s after the bipolar forceps contact with the patient, in order to improve and to favor the process of coagulation.
Visualization of the power that can be activated by foot-switch.

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Two output channels, reciprocally locked, with prevalence of the first activated. The connection of two handpieces with switch button and of a double pedal allows the personalized selection of the functions.
Channel CUT/COAG1 can deliver, though commands from handpiece, the cut currents or the coagulation currents available in COAG1.
Channel CUT/COAG2 can deliver, with commands by handpiece or by foot-switch, the cut currents or the coagulation currents available in COAG2.
The bright signal OVT points out exceeding of the maximum time, recommended for delivery of current, to avoid excessive heating of the treated area.

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Output for the connection of the bipolar accessories with command of activation through double pedal for the selection of the functions.
The bright signal OVT points out exceeding of the maximum time, recommended for delivery of current, to avoid excessive heating of the treated area.

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WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.

PINZE BIPOLARI ~ BIPOLAR FORCEPS ~ PINCES BIPOLAIRES ~ PINZAS BIPOLARES ~ BIPOLÁRNÍ PINZETY
AUTOCLAVABLE 134°

8310110
DIRITTE • STRAIGHT • DROITE • DERECHO • ROVNÁ
11,5 cm

8310115
CURVA • ANGLED • COUDÉE • ANGULO • ZAHNUTÁ
11,5 cm

8310120
DIRITTE • STRAIGHT • DROITE • DERECHO • ROVNÁ
15 cm


8310140
DIRITTE • STRAIGHT • DROITE • DERECHO • ROVNÁ
20 cm

8310122
CURVA • ANGLED • COUDÉE • ANGULO • ZAHNUTÁ
15 cm

8310142
CURVA • ANGLED • COUDÉE • ANGULO • ZAHNUTÁ
20 cm

8310170
CURVA • ANGLED • COUDÉE • ANGULO • ZAHNUTÁ
18 cm

8310175
CURVA • ANGLED • COUDÉE • ANGULO • ZAHNUTÁ
18 cm

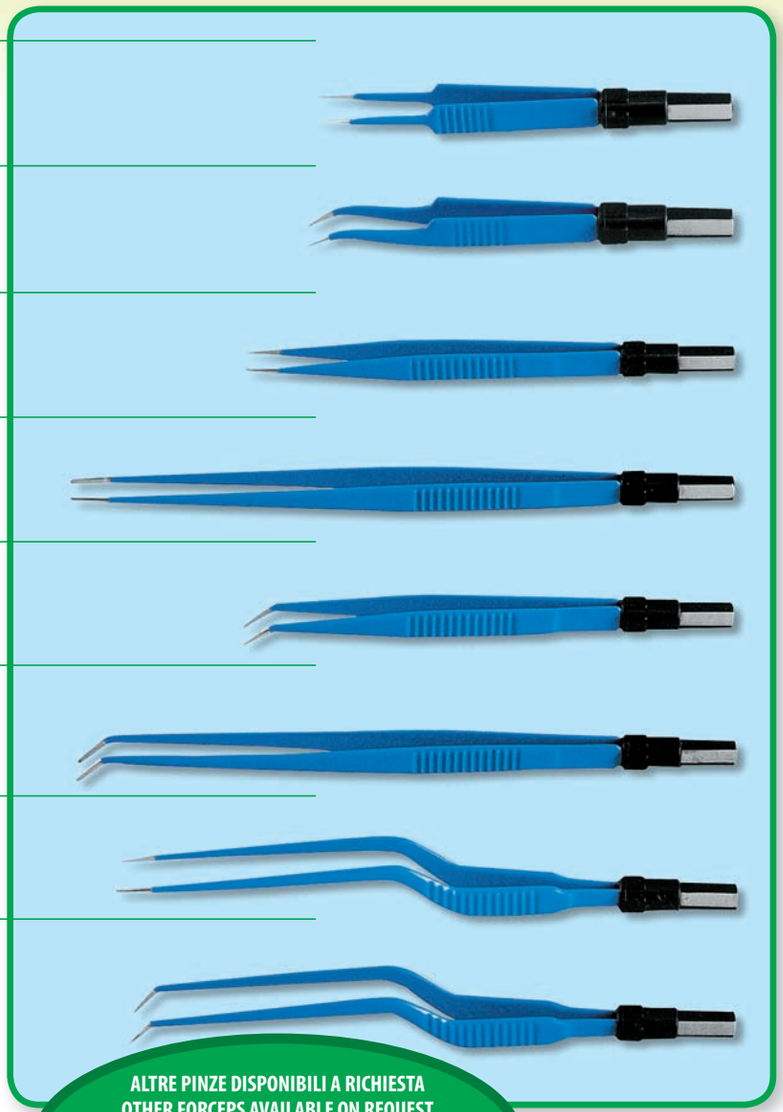ALTRE PINZE DISPONIBILI A RICHIESTA
OTHER FORCEPS AVAILABLE ON REQUEST
AUTRES PINCES DISPONIBLES
OTRAS PINZAS DISPONIBLES
JINÉ DRUHY PINZET K DISPOZICI NA POŽÁDÁNÍ

CE 0120

electrosurgery


8190120
Cavo in silicone per pinze bipolari lungo 3 mt. con connettore da 4 mm. autoclavabile a 134° C.
Silicone cable autoclavable 134° for bipolar forceps 3 m lenght, 4 mm pin
Cordon autoclavable à 134° pour pinces bipolaires.long.3 mt. adaptateur 4 mm.
Cable autoclavable a 134º para pinzas bipolares. largo 3 m., conector 4 mm.
Silikonový kabel autoklávovatelný při 134 °C, pro bipolární pinzety, délka 3 m, 4 mm konektor kolík

00498.04
Adattatore per pinze bipolari
Adapter for bipolar forceps
Adaptateur pour pinces bipolaires
Adaptador para pinzas bipolares.
Adaptér pro bipolární pinzety

WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.

ELETTRODI ED ACCESSORI PER ELETTROBISTURI
ELECTRODES AND ACCESSORIES FOR ELECTROSURGERY UNIT
ELECTRODES ET ACCESSOIRES POUR ELECTROCHIRURGIE
ELECTRODOS Y ACCESORIOS PARA ELECTROCIRUGIA
ELEKTRODY PRO ELEKTROCHIRURGICKÉ JEDNOTKY

CE

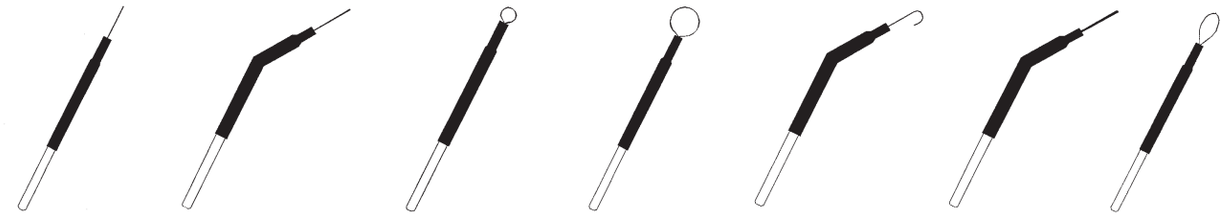00500.L1
00500.L1L

00500.L2
00500.L2L

00500.L3
00500.L3L

00500.L4
00500.L4L

00500.L5
00500.L5L

00500.L6
00500.L6L

00500.L7
00500.L7L

00500.L8
00500.L8L

00500.L9
00500.L9L

00500.L10
00500.L10L

00500.L11

0210
0210/L

0220

0230
0230/L

ELETTRODI • ELECTRODES • ÉLECTRODES • ELECTRODOS ELEKTRODY, DÉLKA 5 cm	
500500.L1	fine diritto • fine straight • très fine droite • fino y recto • rovná
500500.L2	fine curvo • fine angled • très fine coudée • fino angulado • zahnutá
500500.L3	chiuso • loop • lacet • lazo • smyčková, prům **4** mm Ø
500500.L4	chiuso • loop • lacet • lazo • smyčková, prům **8** mm Ø
500500.L5	uncinato curvo • hook angled • crochet coudée • gancho angulado • zahnutá s háčkem
500500.L6	grosso curvo • thick angled • grosse coudée • angulo angulado • silná, zahnutá
500500.L7	chiuso • loop • lacet • lazo • smyčková
500500.L8	chiuso • loop • lacet • lazo • smyčková
500500.L9	sferico • ball • boule • pala • balónková, prům. **3** mm Ø
500500.L10	sferico • angled ball • boule coudée • angulo y con bola • balónková zahnutá prům. **3** mm Ø
500500.L11	microchirurgia • microsurgery • microchirurgie • microcirugía • mikrochirurgická
0210	lancia • lancet • lancette • lanceta • skalpelová
0220	punta • needle • aiguille • aguja • jehlová
0230	sferico • ball • boule • bola • balónková

ELETTRODI • ELECTRODES • ÉLECTRODES • ELECTRODOS ELEKTRODY, DÉLKA 10 cm	
500500.L1L	fine diritto • fine straight • très fine droite • fino y recto • rovná
500500.L2L	fine curvo • fine angled • très fine coudée • fino angulado • zahnutá
500500.L3L	chiuso • loop • lacet • lazo • smyčková, prům **4** mm Ø
500500.L4L	chiuso • loop • lacet • lazo • smyčková, prům **8** mm Ø
500500.L5L	uncinato curvo • hook angled • crochet coudée • gancho angulado • zahnutá s háčkem
500500.L6L	grosso curvo • thick angled • grosse coudée • angulo angulado • silná, zahnutá
500500.L7L	chiuso • loop • lacet • lazo • smyčková
500500.L8L	chiuso • loop • lacet • lazo • smyčková
500500.L9L	sferico • ball • boule • pala • balónková, prům. **3** mm Ø
500500.L10L	sferico • angled ball • boule coudée • angulo y con bola • balónková zahnutá prům. **3** mm Ø
500500.L11L	microchirurgia • microsurgery • microchirurgie • microcirugía • mikrochirurgická
0210L	lancia • lancet • lancette • lanceta • skalpelová
0220L	punta • needle • aiguille • aguja • jehlová
0230L	sferico • ball • boule • bola • balónková


electrosurgery

00500.L0

- Serie completa di 10 pezzi corti
- Complete set short 10 pieces
- Jeu complet court 10 pièces
- Set surtido corto de 10 piezas
- Sada 10 ks krátkých elektrod

00500.L0L

- Serie completa di 10 pezzi lunghi
- Complete set long 10 pieces
- Jeu complet long 10 pièces
- Set surtido largo de 10 piezas
- Sada 10 ks dlouhých elektrod

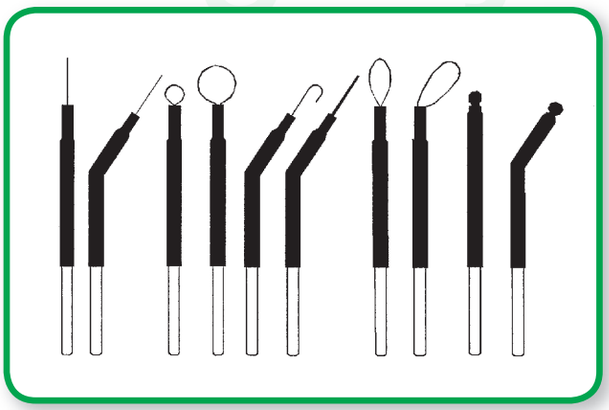WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.

MEDI-PENCIL

- Manipolo portaelettrodi con doppio interruttore cut/coag, non sterile, autoclavabile.
- Electrosurgery pencil with double finger switch, non sterile, autoclavable
- Manche porte electrode réutilisable avec interrupteur digital, non sterile, autoclavable
- Mango porta electrodo interruptor digital, no esteril, esterilizable
- Elektrochirurgický držák monopolárních nástrojů se dvěma spínači, nesterilní, autoklávovatelný

STANDARD TYPE CUT-COAG *surgery pencil*	
2835P300LA	3 metre cable, with blade electrode
2835P300ST	3 metre cable, without electrode
2835P500LA	5 metre cable, with blade electrode
2835P500ST	5 metre cable, without electrode


Spare electrodes for CUT-COAG or foot-controlled electrosurgery pencil	
2835RI00LA	stainless-steel insulated blade electrode
2835RI00AG	stainless-steel insulated needle electrode
2835RI00PA	stainless-steel insulated ball electrode
	(Pack of 10 pcs)

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
DISPOSABLE ELECTROSURGERY PENCIL with CUT-COAG	
28311300LA	electrode and 3 metre cable

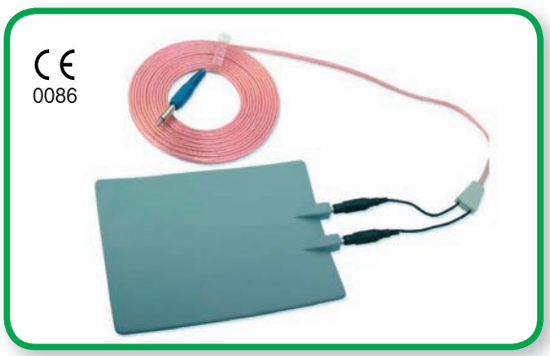PIASTRA IN GOMMA ELETTROCONDUTTIVA
NEUTRAL ELECTRODE OF CONDUCTIVE RUBBER
ELECTRODES NEUTRES EN CAOUTCHOUC CONDUCTEUR
ELECTRODOS NEUTROS DE CAUCHO CONDUCTIVO
NEUTRÁLNÍ ELEKTRODA Z VODIVÉ GUMY

AUTOCLAVABLE 134°

2818C5001A	with bipolar 6,3 mm jack plug
2818C500NP	without connector

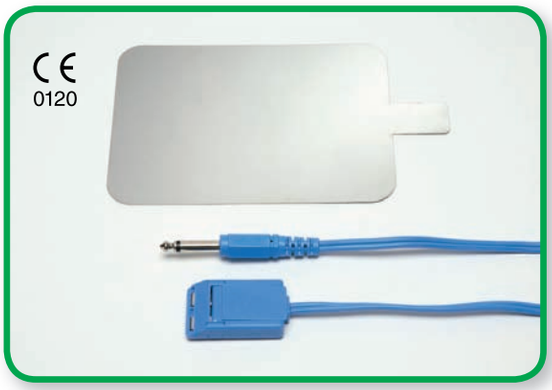00404.08

- **Cavo per piastra**
- **Cable for plate**
- **Câble pour plaque**
- **Cable para placa**
- **Kabel pro neutrální elektrodu**

5365

- **Piastra di riferimento**
- **Neutral plate**
- **Plaque de référence**
- **Placa neutra**
- **Neutrální plochá elektroda**

WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.

MISURE • SIZE • DIMENSIONS • MEDIDAS • ROZMĚRY	
02.860.10	**60** x **40** x **80** H.cm
02.861.10	**70** x **50** x **80** H.cm
02.862.10	**90** x **60** x **80** H.cm

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UTILITY TROLLEY

Aluminium alloy frame. Two removable shelves with upstands on four sides. Complete with one drawer. Mounted on casters Ø 80 mm.
• Knock down construction.

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02.575.00 ~ ELECTROSURGICAL UNITS TROLLEY
Constructed from epoxy-painted steel tube.

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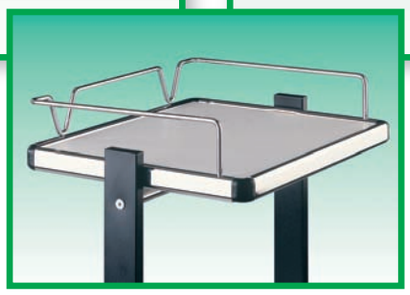
INSTRUMENT TROLLEY

Knock down frame epoxy painted, mounted on casters.
Removable shelves with upstands on three sides.

02.512.00
2 shelves 46x37x77 h cm.
02.513.00
3 shelves 46x37x85 h cm.
02.514.00
4 shelves 46x37x112 h cm. Complete with electrical plug.

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WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.

EC CERTIFICATE

Certificate No 949/MDD

Full Quality Assurance System Approval Certificate

On the basis of our examination carried out according to Annex II, section 3 of Legislative Decree of 1997-02-24, No 46, national transposition of the Directive 93/42/EEC, we hereby certify that:

FAZZINI SRL

manages in the factories of:

VIMODRONE (MI) - S.S. PADANA SUPERIORE 317 (ITA) - Italy

a full quality assurance system ensuring the conformity of the following products:

High frequency surgical equipment

Type ref. DIATROM 50D; DIATROM 80; DIATROM 80D; DIATROM 80F; DIATROM 120; DIATROM 120F; DIATROM 160; DIATROM 160F; DIATROM 200; DIATROM 200F, DIATROM 240; DIATROM 380.

with the relevant essential requirements of the aforementioned national legislation transposing the Directive 93/42/EEC, from design to final inspection and testing.

Reference to IMQ files Nos: 10AG00122.

This Approval Certificate is issued by IMQ S.p.A. as Notified Body for the Directive 93/42/EEC.
Notified Body notified to European Commission under number: 0051.

Date: 2006-07-24

IMQ

This Approval Certificate is subjected to the provisions laid down in the "Rules for managing the EC Certification of Medical Devices on the basis of the Directive 93/42/EEC".
In any case, it does not remain valid after 2011-07-23 (article 11, clause 11 of the Directive).

This is a translation of the Italian text, which prevails in case of doubts

mod.6873/2006-05/1.500/Med.

IMQ S.p.A. I-20138 Milano - Via Quintiliano 43 - tel. 0250731(r.a.) - fax 0250991500 - info@imq.it - www. imq.it
Rea MI 1595884 - Registro Imprese MI 12898410159 - C.F./P.I.: 12898410159 - Capitale sociale 4.000.000 euro.

WE RESERVE THE RIGHT TO CHANGE TECHNICAL SPECIFICATIONS WITHOUT PRIOR NOTICE.